





NATIONAL LIBRARY OF MEDICINE  
Bethesda, Maryland




033173



# TO REDUCE A DISLOCATION OF THE SHOULDER.

Buckle short belt (Q) around the arm above the elbow, one loop being placed on the upper & another on the under side of the arm; through loops pass cords (H) button roll (H) tensely across the shoulder fork (Fig 5); insert (K) of the fork into the socket (D) of counter-bar; apply roll (H) in arm pit; tie cords (H) to foot of the rack (A) strap (M) around the end of roll (H); pass over shoulder and around the opposite end & buckle it tight; apply hollow part of strap (N) to acromion process; pass around under opposite arm and buckle it tight; apply lever and make gradual extension with free motion of arm. NB see pages 62 71 & 89 of Jarvis' Lectures.

THE ADJUSTER.  
MANUFACTURED BY  
GEO. KELLOGG,  
BIRMINGHAM, DERBY, CONN.

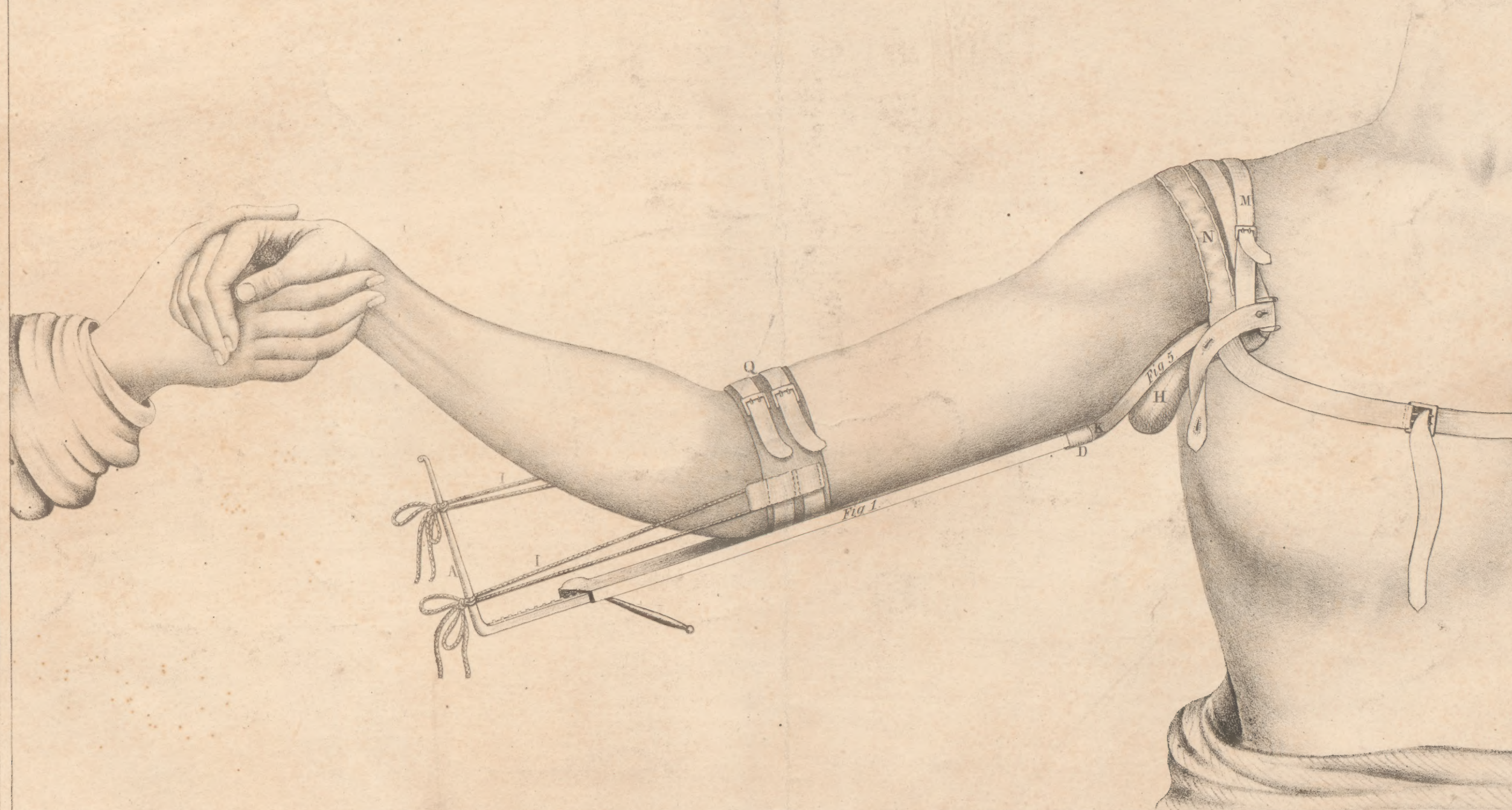

LITHO BY KELLOGG & COMSTOCK, 186 MAIN ST. HARTFORD, CONN.


29858.



# TO DRESS A FRACTURE OF THE LEG.

Place the pads (Fig 5) on each side of limb below the knee; buckle long belt (P) around limb & pads with a loop lying on each side and passing over the head of the pads; arm the loops with cords; (U) apply short belt (Q) around the ankle (or use fracture shoe if you desire it); fasten cords to foot of the rack (A) take out one blade of the jointed fork (Fig 6) insert (K) into its socket (b) of counter-bar; tie cords; (L) to the blade; prepare the plane - place the limb upon it - make the extension required - place bolster (R R) between instrument & plane - place splints on each side of limb, and secure the whole by two tapes (T) or any other bandage. N.B. See page 92 of Larri's Lectures.

THE ADJUSTER,  
MANUFACTURED BY  
GEO. KELLOGG,  
BIRMINGHAM, (DERBY) CONN.

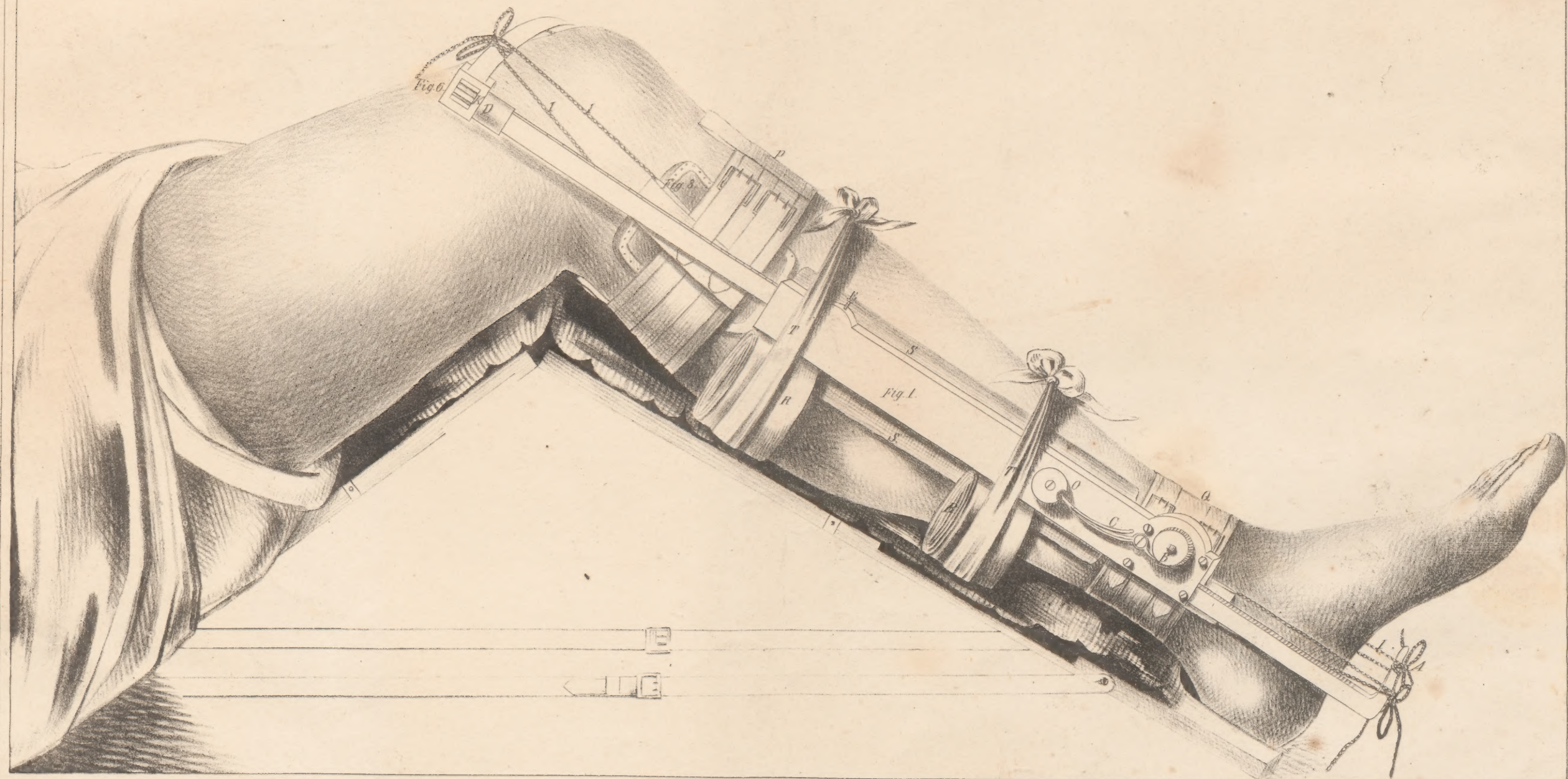


LITH. OF KELLOGG & COMPANY, 130 MAIN ST. HARTFORD, CONN.



# TO DRESS FRACTURES OF THE NECK OF THE THIGH BONE.

Apply long belt (P) without the Pads, roll two straight splints 2 1/2 inches broad, in a square piece of cotton cloth, so as to make a trough, in which place the limb; run out the counterbar (D) and fasten it by screw (E) also the rack bar, so that the foot (A) may cross the sole of the patients foot; use Thigh fork (Fig 4) Apply Instrument arranged on the outside of the inside splint as below; tie cords (H) to foot of the rack; make the necessary extension. Apply hip strap (M) around both arms of the fork; between strap & body, place a thick bolster, so that by tightening strap (M) you may give lateral extension to the neck of the bone. N.B. see page 52 & 92 of Jarvis' Lectures.

THE ADJUSTER,  
MANUFACTURED BY  
GEO. KELLOGG,  
BIRMINGHAM, (DERBY) CONN.

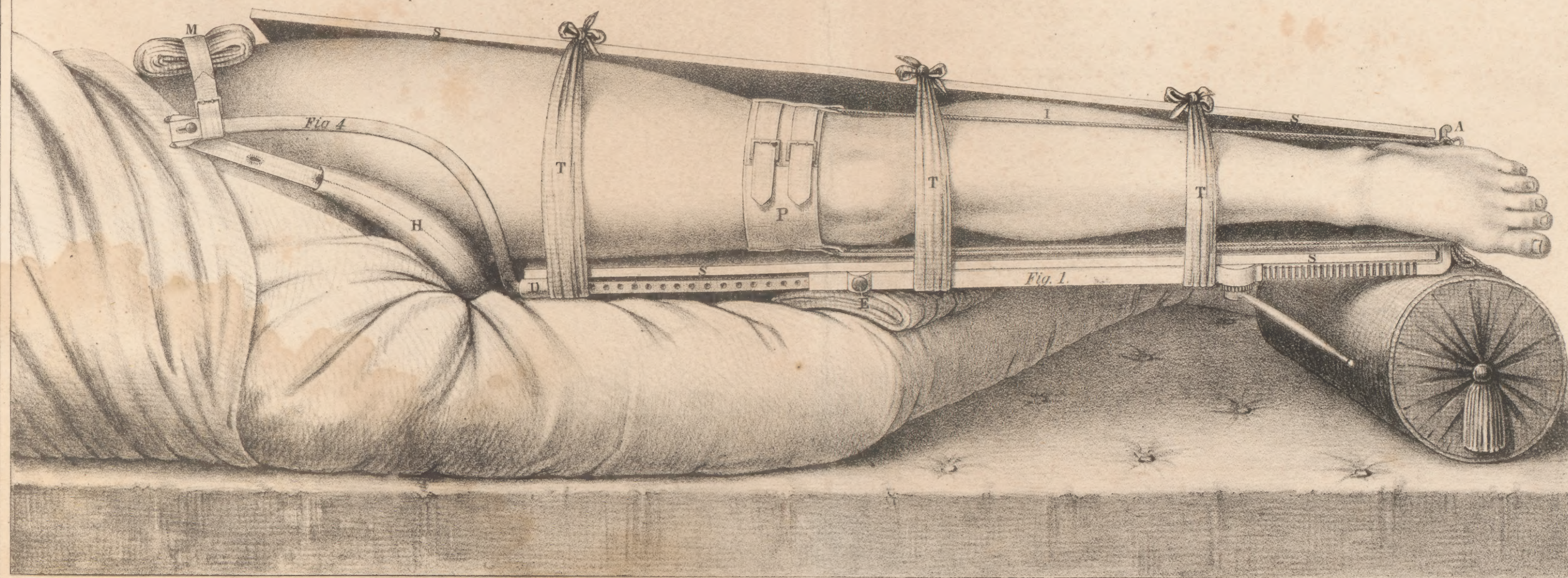


11TH OF KELLOGG & LOWSTOCK 136 MAIN ST. HARTFORD CONN.



# TO DRESS A FRACTURE OF THE LEG.

Place the pads (fig 5) on each side of limb below the knee; buckle long belt (P) around limb & pads with a loop tying on each side and passing over the head of the pads; arm the loops with cords (U) apply short belt (Q) around the ankle (or use fracture shoe if you desire it) fasten cords to foot of the rack (S) take out one blade of the jointed fork (fig 6) insert K into its socket (D) of counter-bar; tie cords (U) to the blade; prepare the plane—place the limb upon it—make the extension required—place bolster (R R) between instrument & plane—place splints on each side of limb, and secure the whole by two tapes (T) or any other bandage. N.B. see page 92 of Jarvis' Lectures.

THE ADJUSTER,  
MANUFACTURED BY  
C. E. KELLOGG,  
BIRMINGHAM, (DERBY) CONN.

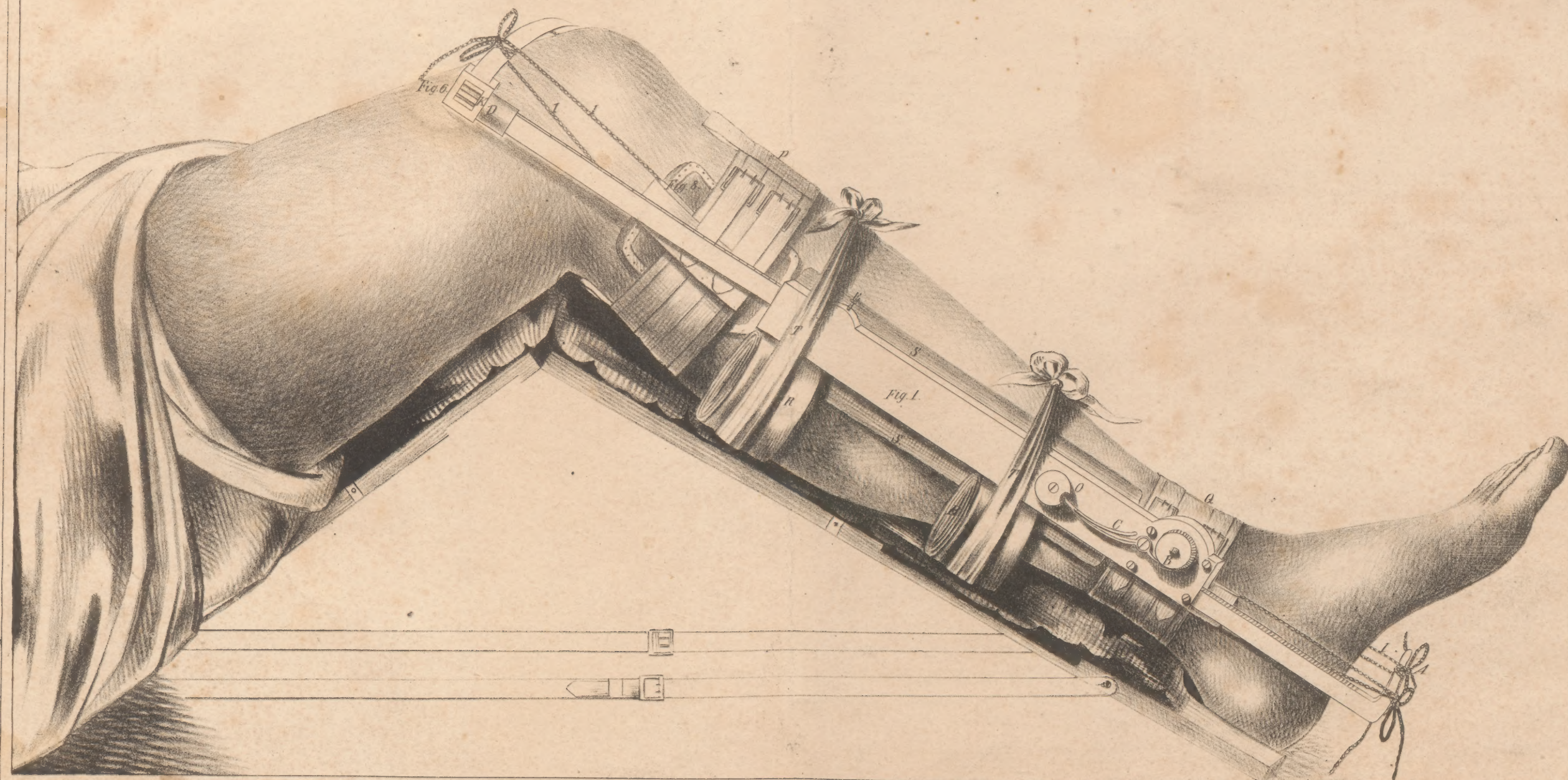




# TO DRESS FRACTURES OF THE NECK OF THE THIGH BONE.

Apply long belt (P) without the Pads, roll two straight splints 2 1/2 inches broad, in a square piece of cotton cloth, so as to make a trough, in which place the limb; run out the counterbar (D) and fasten it by screw (E) also the rack bar, so that the foot (A) may cross the sole of the patients foot; use Thigh fork (Fig 4) Apply Instrument arranged on the outside of the inside splint as below; tie cords (H) to foot of the rack; make the necessary extension. Apply hip strap (M) around both arms of the fork; between strap & body, place a thick bolster, so that by tightening strap (M) you may give lateral extension to the neck of the bone. N.B. see page 52 & 92 of Jarvis' Lectures.

THE ADJUSTER,  
MANUFACTURED BY  
GEO. KELLOGG,  
BIRMINGHAM, (DERBY) CONN.

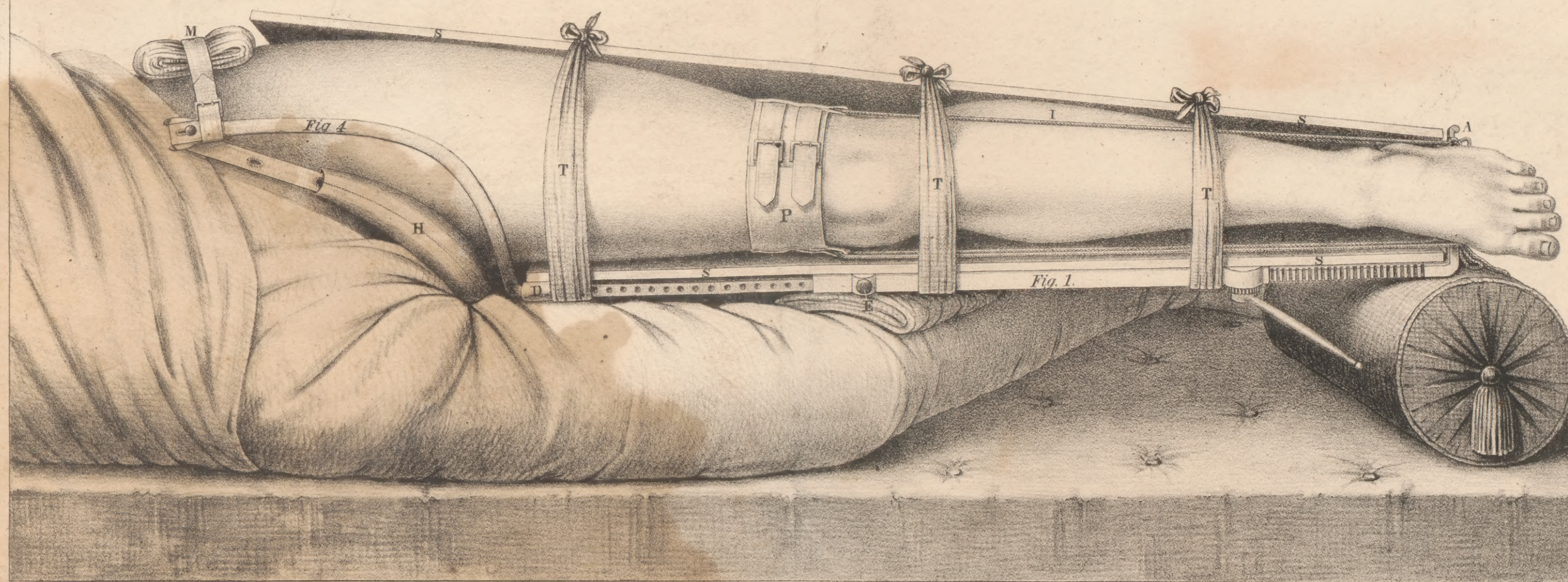

LETTER OF KELLOGG & COMSTOCK, 138 MAIN ST. HARTFORD, CONN.







